# Heated Head Attachment

EN Operating Instructions BG Инструкции за експлоатация  
CS Návod k použití DA Betjeningsvejledning DE Bedienungsanleitung  
EL Οδηγίες λειτουργίας ES Instrucciones de uso ET Kasutusjuhised  
FI Käyttöohjeet FR Mode d'emploi HE הלעפה תוארוה HR Upute za rad  
HU Használati útmutató IT Istruzioni per l'uso JA 取扱説明書  
KO 작동 지침 LT Naudojimo instrukcijos MS Arahan operasi  
NL Gebruiksaanwijzing PL Instrukcje użytkowania  
PT Instruções de operação PT-BR Instruções de operação  
RO Instrucțiuni de utilizare RU Инструкции по эксплуатации  
SK Návod na použitie SL Navodila za uporabo  
SV Driftanvisningar TH คู่มือการใช้งาน ZH 使用说明书  
ZH-TW 繁體中文 操作說明

**EN**

TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, OR PROPERTY DAMAGE, THIS DEVICE MUST BE USED IN ACCORDANCE WITH THE FOLLOWING WARNINGS, CAUTIONS, AND SAFETY INSTRUCTIONS.

# IMPORTANT SAFETY INSTRUCTIONS - ORIGINAL INSTRUCTIONS

Read all safety warnings and instructions before using this device. Failure to follow the warnings and instructions may result in electric shock, fire and/or serious injury.

## CAUTION

- Use this device only for its intended use as described in this manual.
- This device should never be left unattended when plugged in. Always unplug this device from the electrical outlet immediately after charging, prior to use and before cleaning. Unplug from outlet when not in use, and before putting on or taking off parts.
- Suitable for dry location **ONLY**. Do not use outdoors or on wet surfaces and do not expose to water or rain.
- Do not modify or attempt to repair the device. Contains no user-serviceable parts. Please contact Hyperice for service.

- This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge of how to use this device, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
- The device has a heated surface. Persons who are sensitive to heat must be careful when using the device.
- Children should be supervised to ensure that they do not play with the device.
- Unplug the charger from battery after charging or prior to use. Heated Head will not function while charger is plugged into battery.
- Never leave the Heated Head operating or charging unattended.
- Do not immerse device in water. Clean only with a damp cloth.
- Do not use Heated Head on wet skin or over topical ointment, especially ones that contain heat producing ingredients as skin burns could result.
- Do not attempt to disassemble.

#### INTERNAL BATTERY INFORMATION

The Heated Head is equipped with a rechargeable lithium polymer battery. The internal battery is designed to allow use of the Heated Head anywhere—even when power outlets aren't available. The Heated Head may need to be charged before first use. A fully charged battery will provide power for 60+ minutes of continuous use. It takes approximately 3.5 hours to fully charge the battery when the device is plugged in and not in use.

## WARNING

- Use heated surfaces carefully. May cause serious burns. Do not use over insensitve skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitve skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the device frequently to reduce the risk of blistering. Immediately remove from body and turn off at first sign of discomfort.

- Do not charge the Heated Head overnight or leave the Heated Head unattended when charging or in use.
- This device contains a Lithium-polymer battery. Batteries are only replaceable by a skilled person with a special tool. Contact Hyperice for service or repair. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of. The battery must be complied with safely at an appropriate e-waste disposal or recycling facility.
- Follow all charging instructions and do not charge the battery pack or device outside of the temperature range specified in the instructions. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
- Nominal temperature range for charging the battery is 0°C ~ 45°C (32°F ~ 113°F).
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. If contact accidentally occurs, flush with water. If liquid contacts eyes, seek medical help. Liquid ejected from the battery may cause irritation or burns.
- The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, do not expose a battery pack or device to fire or excessive temperature. Exposure to fire or temperature above 130°C (265°F) may cause explosion.
- In the event of the battery pack leaking, do not allow any liquid to come into contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical attention.
- In case of battery rupture, major leakage or explosion, evacuate the area and quarantine the contaminated area. Provide good ventilation to clear out any fumes, gases or pungent odors.
- Eyes – Rinse eyes with plenty of water for 15 minutes; Seek immediate medical attention.
- Skin – Rinse affected area with plenty of water and soap or take a shower for 15 minutes.
- Inhalation – Expose the person to fresh air and use artificial respiration if needed; Seek medical attention if necessary.
- Ingestion – Consult a physician or local poison control center immediately
- Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

- Avoid applying the heated head to the same location for a prolonged period of time. Exposing skin to heated surfaces for a long period of time may cause discomfort or burns.

#### CHARGING


To charge, connect the USB-C cord into the battery's charging port and plug the cord into a USB-A adapter. Then plug the adapter into the wall. Fully charge the battery for up to three hours before first use. A LED light pulses to show the battery level and indicate active charging. Full charge is indicated when the LED light stops pulsing and remains fully illuminated green. Pressing the power button will turn on the device and will display the battery level. Device may be recharged at any time and at any battery level. It is not recommended to fully deplete the battery to the orange LED level. Recharge the device monthly to best prolong battery life.

#### OPERATING

Ensure battery is charged. On the heated head, press and hold the power button for 0.5s to activate the device. The battery led will illuminate to show the Heated Head is activated. A single press on the power button will then cycle between level 1, 2, and 3 heat levels. You can adjust the heat levels at any time. The device will turn off automatically after 15 minutes without further input. To turn off the device, press and hold the power button for 0.5s. The average run-time is around 60 minutes at level three heat, and longer at lower levels.

#### MAINTENANCE, CLEANING, AND STORAGE

Use a towel or soft cloth to wipe any moisture from device. Keep device away from heat sources, liquids, children, and pets.
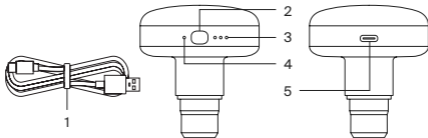
## SPECIFICATIONS

1. USB-C cord
2. Power button (ON/OFF) / heat setting button
3. Heat level indicators
4. Battery life indicator
5. USB-C charging port

Heat: Level 1—109°F (43°C), Level 2—115°F (46°C), Level 3—120°F (49°C)

Weight: Heated Head—3.2 oz., USB-C cord—0.5 oz.

Electrical requirement: ⚡ 5V ===, 1A

## RECYCLING

Please visit <https://hyperice.com/recycling> for information about recycling this packaging.

# CAUTION

Not intended for medical use in hospitals. Intended for household or commercial use. When not in use, it is recommended to store the device between -10°C and 25°C.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

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#### WARRANTY

This product is covered by a limited warranty from Hyperice. Please visit [hyperice.com/warranty](https://hyperice.com/warranty) to review the warranty in your country.

## FOR MEXICO ONLY

- This device is not intended for use by persons (including children) whose physical, sensory or mental capabilities are different or reduced, or who lack experience or knowledge, unless such persons have been given supervision or training in the operation of the device by a person responsible for their safety.
- Children should be supervised to ensure that they do not use the devices as toys.
- The device must only be powered at the safety extra-low voltage indicated on the device.
- The battery must be removed from the device before disposal. The device must be disconnected from the charger when the battery is removed. The battery is removed for safe disposal.
- Operation of this equipment is subject to the following two conditions:  
(1) this equipment or device may not cause harmful interference, and  
(2) this equipment or device must accept any interference, including interference that may cause undesired operation.

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Copyright 2023 © Hyperice, Inc. All rights reserved.
Hyperice, Inc., 525 Technology Drive, Suite 100