Normatec Lower Legs

EN Operating Instructions AR لي غش التا ادي لعت
BG Инструкции за експлоатация CS Návod k použití DA Betjeningsvejledning
DE Bedienungsanleitung EL Οδηγίες λειτουργίας ES Instrucciones de uso
ET Kasutusjuhised FI Käyttöohjeet FR Mode d'emploi HE הלעפע תוארה
HR Upute za rad HU Használati útmutató IT Istruzioni per l'uso
LT Naudojimo instrukcijos MS Arahan operasi NL Gebruiksaanwijzing
PL Instrukcje użytkowania PT Instruções de operação
PT-BR Instruções operacionais RO Instrucțiuni de utilizare SK Návod na použitie
SL Navodila za uporabo SV Driftanvisningar

EN

TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE, AND PERSONAL INJURY, OR PROPERTY DAMAGE, THIS DEVICE MUST BE USED IN ACCORDANCE WITH THE FOLLOWING WARNINGS, CAUTIONS, AND SAFETY INSTRUCTIONS.

IMPORTANT SAFETY INSTRUCTIONS - ORIGINAL INSTRUCTIONS

Read all instructions before using the Normatec Lower Legs.

DANGER

To reduce the risk of electric shock:

- Always unplug this device from the electrical outlet immediately after using and before cleaning.
- Do not reach for a device that has fallen into water. Unplug immediately.
- Do not use while bathing or in a shower.
- Do not place or store device where it can fall or be pulled into a tub or sink. Do not place in or drop into water or other liquid.

WARNING

To reduce the risk of burns, fire, electric shock, or injury to persons:

- A device should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- This device can be used by children aged from 8 years and above and persons with reduced physical, sensory, or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children should be supervised to ensure that they do not play with the appliance.

- Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- Never operate this device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or if it has gotten wet. Return the device to a service center for examination and repair.
- Do not carry this device by supply cord or use cord as a handle.
- Keep the cord away from heated surfaces.
- Do not use outdoors.
- Do not operate where aerosol (spray) devices are being used or where oxygen is being administered.
- Cleaning and user maintenance shall not be done by children without supervision.
- Do not remove screws or attempt to disassemble.
- Unplug the device after charging and prior to use.
- Use device only according to Operating Instructions.
- When using the Normatec Lower Legs use only on dry, clean surfaces of the body. Please immediately stop use of the Normatec Lower Legs upon experiencing any pain or discomfort. Use of the Normatec Lower Legs should also be avoided on all bruises, contusions, rashes or irritated or injured areas of the skin until they are fully healed.

 This symbol indicates this device is intended to be supplied from a detachable supply unit for purposes of recharging the battery. The Normatec Lower Legs is intended to be charged by a 5V USB. Do not charge the Normatec Lower Legs overnight or leave the Normatec Lower Legs unattended when charging or in use.

Please do not use the Normatec Lower Legs without first obtaining approval from your physician if any of the following apply:

- Pregnancy, diabetes with complications such as neuropathy or retinal damage, wear of pacemakers, recent surgery or injury, epilepsy or migraines, herniated disks, spondylolisthesis, spondylolysis, or spondylosis, recent joint replacements or IUD, metal pins or plates or any concerns about your physical health. Frail individuals and children should be accompanied by an adult when using any percussion or vibration device. These contraindications do not mean that you are not able to use a percussion or vibration device but we advise that before doing so you consult your physician. Ongoing research is being done on the effects of percussion massage for specific medical disorders that may result in the shortening of the list of contraindications as shown above.

The Normatec Lower Legs contains a battery which must be complied with safely at an appropriate e-waste disposal or recycling facility. This device contains batteries that are non-replaceable. This appliance contains batteries that are only replaceable by skilled persons.

CAUTION

The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not disassemble, heat above 100° C or incinerate. Dispose of used battery promptly. Keep away from children. Do not disassemble and do not dispose of in fire.

SAVE THESE INSTRUCTIONS

The Normatec Lower Legs is an air massage device that delivers targeted pressure to the calves to care for muscles, relieve tension, provide a relaxing massage, accelerate warmup and recovery, and help maintain flexibility and range of motion.

ILLUSTRATIONS

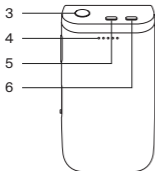
Normatec Lower Legs Device (single-person use only)

1. HyperSync™ indicator
2. Bluetooth® connection indicator
3. Start / Stop button
4. Battery level indicator
5. Time adjustment button
6. Pressure level adjustment button
7. Power (ON/OFF) button
8. USB-C charging port
9. Leg compression wrap
10. USB-C cord (2 included)


TOP



FRONT

SIDE

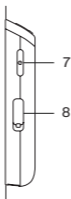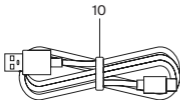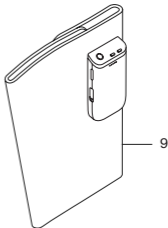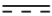
FULL DEVICE

LABELS

The following labels and symbols appear on the control unit, attachments, and/or packaging.

| Symbol | Description | Location |
|-------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|---------------------------|
| IP22 | Degree of protection against ingress of water | On fabric tag |
|  | Direct current | On charger |
|  | Alternating current | In manual |
|  | Manufacturer's name and address | On fabric tag |
|  | Separate collection for waste electrical and electronic equipment. | On fabric tag |
|  | Serial number of the console | On base of control unit |
|  | Keep dry | On package and fabric tag |
|  | Place in and out of standby mode | On side of control unit |
|  | Do not wash | On fabric tag |
|  | The Bluetooth figure mark | On packaging |
|  | Warning symbol to identify a hazard that may lead to death or serious injury | In manual |
|  | Caution symbol to indicate the need for the user to consult the instructions for use for important cautionary information such as warnings and precautions that cannot, for a variety of reasons, be presented on the medical control unit itself | In manual |
|  | Tip to provide guidance to make use easier. Risk to user is considered negligible | In manual |
|  | Catalog/Model Number | On fabric tag |
|  | Made in China | On fabric tag |
|  | Recycle Lithium-Ion | On control unit |

OPERATING INSTRUCTIONS

CHARGING THE SYSTEM

- Step 1: Plug the USB cable into the control unit
- Step 2: Plug the USB cable into an USB-A 5V source
- While the control unit is charging the battery status indicators will cycle from left to right. When all 5 battery status indicators are on ON, this means that the system is fully charged.

The control unit will NOT power on and cannot be operated while the USB cable is connected.

Device is intended to be charged with the power adaptor provided or by 5V USB. If you have any questions, please contact Hyperice at +1.949.565.4994 or your local Customer Care Center or representative.

SET UP THE SYSTEM

- Step 1: Secure the wrap tightly around your calf. Affix the two sides of the hook and loop material together for a snug but comfortable fit.
- Step 2: Press the power button on the Normatec Lower Legs control unit firmly for one second to turn on the system. While the control unit is on, the battery LEDs and display screen will light up.

USING HYPERSYNC™

Your Normatec Lower Legs devices arrive paired, which means they will inflate and deflate on the same cycle. To unpair devices, press and hold the Start / Stop button on either device until the device screen reads "Unpaired!" and the white pairing indicator light is no longer illuminated. To pair your devices, turn on both devices by pressing and holding the power (ON/OFF) button. Then press and hold Start / Stop button on either control unit until the display screen says "Pairing!" Finally, press and hold the Start / Stop button on the other device until both devices say "Pairing!" and the white pairing indicator light is illuminated.

CONNECT TO THE HYPERICE APP

Download the Hyperice App from the App Store or the Google Play Store. To connect your device to the Hyperice App via Bluetooth® make sure that your device is turned on, Bluetooth® is turned on in your phone, and your device is within range of your phone. The Bluetooth LED on the pod will pulse blue, which means it is waiting to pair. If prompted, tap "Scan for Devices." Select your device when it pops up on the screen. Once you pair your Normatec Lower Legs, HyperSmart™ will allow you to automatically start your session and control the pressure level and time from within the Hyperice App.

CYBERSECURITY

It is recommended to configure the Hyperice App for automatic updates to ensure cybersecurity. It is also recommended to keep your Operating System up to date and to configure your Operating System for automatic updates.

ADJUST THE PRESSURE LEVEL

Adjust the pressure level of the session by pressing the pressure level adjustment button on the top left of the control panel. Pressure level 1 is the gentlest setting. The massage becomes more intense as the pressure level is increased. Level can be adjusted while the session is running.

ADJUST THE SESSION TIME

Adjust the session time by pressing the time adjustment button. The session time can be set to 15, 30, 45, and 60 minutes. Time can be adjusted while the session is running. Press the time adjustment buttons to cycle through the session time options to add or subtract time from the session.

START THE SESSION

To start the session, press the Start/Stop button.

STOP OR PAUSE THE SESSION

To stop the session at any time, press the start/stop button. This will pause the session. To restart your paused session, press the start/stop button again. If you are done using the system, remove the wraps from your limbs, and turn off the control unit by pressing the power button.

FINISH THE SESSION

The session will continue massaging until time runs out and the display reads "Finishing Cycle". The system will continue until the current cycle is finished. When the session is complete, remove the wraps from your limbs, and turn off the control unit by pressing the power button.

TURN OFF THE CONTROL UNIT

To turn off the system, press the power button and confirm that the Battery Status Indicators and Display are off.

CLEANING THE SYSTEM

- Wipe down the inside and out with a damp, clean cloth.
- Dry thoroughly with a clean cloth.
- Do not machine wash or dry.
- Do not dry clean.

MAINTAINING THE SYSTEM

The system requires no routine maintenance or service except for the care in this section.

STORING THE SYSTEM

Store the system in a clean, dry location.

RECYCLING

Please visit <https://hyperice.com/recycling> for information about recycling this packaging.

TECHNICAL INFORMATION

Do not attempt to take apart the system. The system has no user-serviceable parts. There are no user-replaceable fuses.

PRODUCT SPECIFICATIONS

- Model: CLV23
- Dimensions (fully open): 25" (width), 1" (depth), 12" (height); [63.5 cm (width), 2.5 cm (depth), 30.5 cm (height)]
- Weight: 1.2 lbs [0.5 kg]

- Electrical requirement: 5V DC 1A
- Maximum Air Pressure: 220 mm Hg
- Internal Battery Rating: 3.6V, 2550mAh/3.65V, 2600mAh
- Ambient Temperature (charging): + 50 °F to 95 °F [+10 °C to 35 °C]
- Temperature (operating): +41° F to 104° F [+5° C to +40° C]
- Temperature (storage): -13° F to +158° F [-25° C to +70° C]
- Relative Humidity (operating): 15% to 93%, non-condensing
- Relative Humidity (storage): -25° C without relative humidity control; +70° C at relative humidity up to 93%, non-condensing

| | HVIN: V1.0 | HVIN: V1.2 |
|----------------------|---------------------------------------|----------------------------------------|
| Frequency | 2402 to 2480 MHz, 433.92 MHz | 2402 to 2480 MHz, 433.92 MHz |
| Modulations | GFSK, FSK | GFSK, FSK |
| Transmit Power | +3.39 dBm (BLE), 2.05 dBm (ULP) | +3.47 dBm (BLE), 3.71 dBm (ULP) |
| Receiver Sensitivity | -96 dBm (BLE mode), -118 dBm (ULP) | -96 dBm (BLE mode), - 118 dBm (ULP) |
| Security | AES HW | AES HW |

Hereby, NormaTec declares that the radio equipment type CLV23 is in compliance with Directive 2014/53/EU. The full text of the EU declaration of conformity is available at the following internet address: <https://hyperice.com/doc>

TROUBLESHOOTING

| Problem | Possible Causes | Solutions |
|---------------------------|--------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| The system does not start | Power is not turned on, battery is not charged | Press the power button to turn the control unit on. Check that the battery is charged. |
| The wrap does not inflate | The session has not been started The wrap is punctured or damaged | Check that the start / stop LED is ON and the timer is counting down. If not, press start / stop button to start or resume session. Check that there is no air leak error message on the display screen. |

| | | |
|-----------------------------------------------------|----------------------------------------|-------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| The system stopped pumping | The session has been stopped or ended. | <p>Check that the start / stop LED is ON and the timer is counting down. If not, press start / stop button to start or resume session.</p> <p>Check that there is no air leak error message on the display screen.</p> <p>Check that there is no Low Battery message on the display screen.</p> |
| Air leak message: Air Leak | Air leak | Check for punctures in the wrap. |
| Low Battery | Battery needs to be charged | Plug in the control unit to charge the battery. |
| Cannot establish or maintain a Bluetooth connection | Bluetooth is turned off | Turn on Bluetooth on the phone attempting to pair with the Normatec Lower Legs control unit. |
| Devices not syncing settings and/or treatment | Devices are not paired. | Check that the white pairing indicator light is illuminated on both devices. If not, your devices are not paired. Follow the "Using HyperSync™" instructions in the Operating Instructions to pair your devices. |

Call Hyperice customer service at +1.949.565.4994 if further assistance is needed.

WARRANTY

This product is covered by a limited warranty from Hyperice. Please visit hyperice.com/warranty to review the warranty in your country.

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Copyright 2023 © Hyperice, Inc. All rights reserved.